

HRV 3C Protease ELISA Kit

Catalog Number LDG00017E

Package 96 T (8 x 12 strips) / Customized package

For full product information, images and publications, please visit [our website](#).



Overview

Components

Items	Quantity
HRV 3C Protease ELISA plate	96 wells (12 x 8-well strips)
HRV 3C Protease Standard	2 vials (Lyophilized form)
Standard reconstitution buffer	2 vials (1.1 mL)
Standard & Sample diluent buffer	1 vial (12 mL)
HRP-antibody conjugate	1 vial (65 µL)
HRP-antibody conjugated diluent buffer	1 vial (12 mL)
20× wash buffer	1 vial (15 mL)
TMB	1 vial (12 mL)
Stop solution	1 vial (6 mL)
Microplate sealing film	1 sheet

Specifications

Reactivity

Human rhinovirus

Application

Sandwich ELISA analysis

Tainan Headquarter

+886-6-2536677 +886-6-2531536

info@leadgene.com.tw

Taipei Office

+886-2-27065528

Detection Method

Colorimetric

Specificity

HRV 3C Protease

Sample Type

Samples contain HRV 3C Protease

Sensitivity

Limit of detection (LoD): 0.057 ng/mL

Limit of quantification (LoQ): 0.19 ng/mL

Assay Range

1.5625-100 ng/mL

Instruction

Shipping

The product is shipped with polar packs. Upon receipt, store it immediately at 2-8°C for long term storage.


Stability & Storage

This product is stable after storage at:

- 2-8°C for unopened product.

Please refer to product manual for storage constructions.

Image

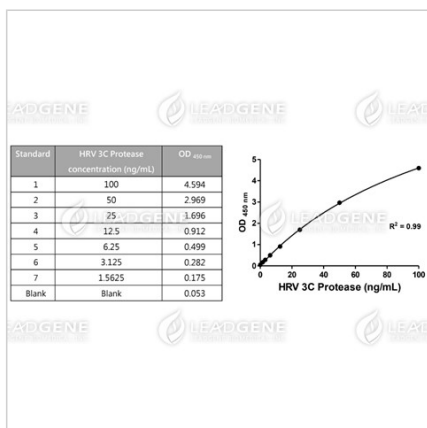

Typical data
The following data are for demonstration only

Reagents (Store at 2-8°C)	Quantity (x96 well kit)	Reconstitution
HRV 3C Protease ELISA plate	96 wells (12 x 8 well strips)	Ready for use
Standard	2 vials (lyophilized form)	Refer to the vial label for reconstitution volume. Reconstitute by adding Standard reconstitution buffer to be a stock solution of 1000 ng/mL. (see procedure, section 8.2)
Standard reconstitution buffer	2 vials (1.1 mL)	Ready for use
Standard & Sample diluent buffer	1 vial (12 mL)	Ready for use
HRP antibody conjugate	1 vial (85 µL)	Dilute 200 x with HRP antibody conjugated diluent buffer (see reagent preparation, section 5.A)
HRP antibody conjugated diluent buffer	1 vial (12 mL)	Ready for use
20 X wash buffer	1 vial (15 mL)	Dilute 20 x with distilled water (see reagent preparation, section 5.B)
TMB	1 vial (12 mL)	Ready for use
Stop solution	1 vial (6 mL)	Ready for use
Microplate sealing film	1 sheet	N/A

Reagents provided and reconstitution

	1	2	3	4
A	Standard 1 (100 ng/mL)	Standard 1 (100 ng/mL)	Sample 1	Sample 5
B	Standard 2 (50 ng/mL)	Standard 2 (50 ng/mL)	Sample 1	Sample 5
C	Standard 3 (25 ng/mL)	Standard 3 (25 ng/mL)	Sample 2	Sample 6
D	Standard 4 (12.5 ng/mL)	Standard 4 (12.5 ng/mL)	Sample 2	Sample 6
E	Standard 5 (6.25 ng/mL)	Standard 5 (6.25 ng/mL)	Sample 3	Sample 7
F	Standard 6 (3.125 ng/mL)	Standard 6 (3.125 ng/mL)	Sample 3	Sample 7
G	Standard 7 (1.5625 ng/mL)	Standard 7 (1.5625 ng/mL)	Sample 4	Sample 8
H	Blank	Blank	Sample 4	Sample 8

An example of orientation of standards, blanks and samples in the stripwells microplate

Disclaimer : For Research Use or Further Manufacturing Only.

Tainan Headquarter

Taipei Office

+886-6-2536677 +886-6-2531536

+886-2-27065528

info@leadgene.com.tw

Tainan Headquarter

☎ +886-6-2536677 ☎ +886-6-2531536

✉ info@leadgene.com.tw

Taipei Office

☎ +886-2-27065528