

LeadGMP® PNGase F

Catalog Number LDG004R-GMP

Package 15 KU / 75 KU / Customized package

For full product information, images and publications, please visit [our website](#).



Overview

Description

PNGase F is an enzyme used in biochemistry and molecular biology to remove N-linked glycans from glycoproteins. By using PNGase F, researchers can enzymatically cleave between these sugar chains and asparagine residues of glycoproteins, allowing for the study of protein structure and function, particularly in glycosylation research.

Components

LeadGMP® PNGase F	1 vial (400 U/μL)
10X Glycoprotein Denaturing Buffer	1 vial (1 mL)
10X Reaction buffer	1 vial (1 mL)
10% NP-40	1 vial (1 mL)

Specifications

Expression system

Escherichia coli

Concentration

400 U /μL

Buffer

20mM Tris-HCl, 50mM NaCl, 5mM EDTA, pH 7.5

Purity

>95% as determined by SDS-PAGE analysis.

Tainan Headquarter

+886-6-2536677 +886-6-2531536

info@leadgene.com.tw

Taipei Office

+886-2-27065528

Unit Definition

One unit is defined as the amount of enzyme required to remove > 95% of the carbohydrate from 10 µg of denatured RNase B in 1 hour at 37°C in a total reaction volume of 10 µl.

Mycoplasma

Not detected

Endotoxin level

<0.05 EU per 1 µg of the protein by the LAL method.

Form

Liquid

Instruction**Stability & Storage**

This product is stable after storage at:

- -20°C for long-term storage under sterile conditions. Avoid repeated free-thaw cycles.

Manufacturing specifications

LeadGMP® recombinant proteins are manufactured in ISO 13485:2016 and GMP certified facility. The processes include:

- Animal-free reagent and laboratory
- Manufactured and tested under GMP guideline
- Testing and traceability of raw material
- Records of the maintenance and equipment calibration
- Personnel training records
- Batch-to-batch consistency
- Documentation of QA control and process changes
- Manufactured and tested under an ISO 13485:2016 certified quality management system
- Stability monitor of product shelf-life


Image**Tainan Headquarter**

+886-6-2536677 +886-6-2531536

info@leadgene.com.tw

Taipei Office


+886-2-27065528

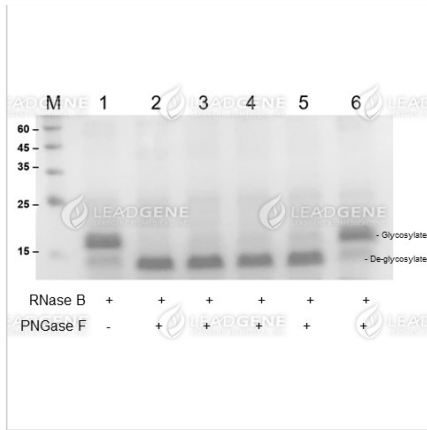

The standard assay was performed by incubating different amount of PNGase F and 10 µg of RNase B.

SDS-PAGE analysis of RNase B digested with LeadGMP® PNGase F .

Lane 1: Uncut, Lane 2: 0.01 U, Lane 3: 0.005 U, Lane 4: 0.0033 U, Lane 5: 0.0025 U, Lane 6: 0.002 U

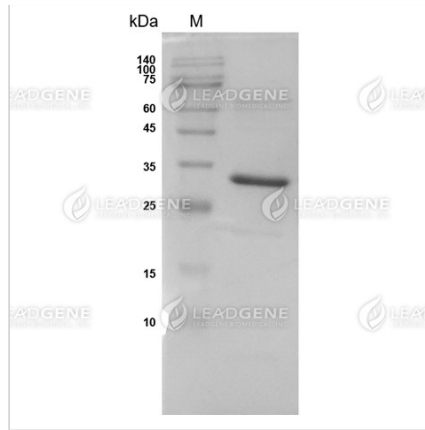

SDS-PAGE analysis of recombinant PNGase F